Database

Open Access

ReplicationDomain: a visualization tool and comparative database for genome-wide replication timing data

Nodin Weddington[†], Alexander Stuy[†], Ichiro Hiratani[†], Tyrone Ryba, Tomoki Yokochi and David M Gilbert*

Address: Department of Biological Sciences, Florida State University, Tallahassee, Florida 32306, USA

Email: Nodin Weddington - weddington@bio.fsu.edu; Alexander Stuy - astuy@bio.fsu.edu; Ichiro Hiratani - ihiratani@bio.fsu.edu; Tyrone Ryba - tryba@bio.fsu.edu; Tomoki Yokochi - yokochit@bio.fsu.edu; David M Gilbert* - gilbert@bio.fsu.edu

* Corresponding author †Equal contributors

Published: 10 December 2008

BMC Bioinformatics 2008, 9:530 doi:10.1186/1471-2105-9-530

This article is available from: http://www.biomedcentral.com/1471-2105/9/530

© 2008 Weddington et al; licensee BioMed Central Ltd.

This is an Open Access article distributed under the terms of the Creative Commons Attribution License (<u>http://creativecommons.org/licenses/by/2.0</u>), which permits unrestricted use, distribution, and reproduction in any medium, provided the original work is properly cited.

Received: I September 2008 Accepted: 10 December 2008

Abstract

Background: Eukaryotic DNA replication is regulated at the level of large chromosomal domains (0.5–5 megabases in mammals) within which replicons are activated relatively synchronously. These domains replicate in a specific temporal order during S-phase and our genome-wide analyses of replication timing have demonstrated that this temporal order of domain replication is a stable property of specific cell types.

Results: We have developed ReplicationDomain <u>http://www.replicationdomain.org</u> as a webbased database for analysis of genome-wide replication timing maps (replication profiles) from various cell lines and species. This database also provides comparative information of transcriptional expression and is configured to display any genome-wide property (for instance, ChIP-Chip or ChIP-Seq data) via an interactive web interface. Our published microarray data sets are publicly available. Users may graphically display these data sets for a selected genomic region and download the data displayed as text files, or alternatively, download complete genome-wide data sets. Furthermore, we have implemented a user registration system that allows registered users to upload their own data sets. Upon uploading, registered users may choose to: (1) view their data sets privately without sharing; (2) share with other registered users; or (3) make their published or "in press" data sets publicly available, which can fulfill journal and funding agencies' requirements for data sharing.

Conclusion: ReplicationDomain is a novel and powerful tool to facilitate the comparative visualization of replication timing in various cell types as well as other genome-wide chromatin features and is considerably faster and more convenient than existing browsers when viewing multi-megabase segments of chromosomes. Furthermore, the data upload function with the option of private viewing or sharing of data sets between registered users should be a valuable resource for the scientific community.

Background

In eukaryotic cells, segments of chromosomes replicate via the synchronous firing of clusters of replication origins

[1]. These segments or "replication domains" replicate in a defined temporal order during S-phase. This replication-timing program is cell type specific [2], and developmen-

tally regulated changes in this program are associated with changes in chromatin structure and gene expression [2-5]. In particular, a global re-organization of this replicationtiming program occurs during the differentiation of mouse embryonic stem cells (mESCs), with changes occurring at the level of large (~600 kb) chromosomal domains reflecting global re-positioning of sequences within the nucleus [2]. Moreover, pluripotent cells can be distinguished from differentiated cells not only by differences in their replication timing profiles but by their smaller and more numerous replication domains [2]. Hence, replication timing is a unique epigenetic property of chromatin in that it is regulated at the level of megabase-sized domains. Establishing replication maps for various tissues is likely to provide a database of chromosome segments that undergo large changes in organization during differentiation.

The significance of a replication-timing program has remained elusive. In several model systems, defects in replication-timing are associated with defects in chromosome condensation, sister chromatid cohesion, and genome stability [6,7]. Abnormal replication-timing control has become a clinical marker for predicting malignant cancers [8-12]. In particular, specific chromosome translocations result in a chromosome-wide delay in replication timing that triggers additional chromosome translocations at a high frequency [13,14]. Cells from patients with several inherited human diseases show defects in replication-timing that correlate with mis-regulation of genes during development [15-18]. Also, replication domains are separated by timing transition regions (the domain boundaries) that appear to be devoid of replication origins, requiring that a single replication fork travel very long distances between early and late replicating domains [2,19,20]. Evidence suggests that genes lying within these transition regions are prone to DNA damage [21,22]. While very few such boundaries have been mapped, their cell-type specificity suggests the possibility that differential organization of replication domains may contribute to cell type specific predispositions to certain types of DNA damage. Hence, establishing a database of replication timing profiles for various tissues and their relationship to transcription and other chromosomal properties is a prerequisite for understanding the roles of replication timing in chromosome-based diseases. These roles may extend beyond epigenetic regulation of transcription: the locations and directions of replication forks, the organization of replication complexes that coordinate replication of large domains, and the locations of domain boundaries may constitute an epigenetic basis for tissuespecific or cancer-promoting differences in genome stability.

Surprisingly few genome-wide studies of replication timing have been performed [23]. Early studies in *Drosophila* cells with cDNA arrays [24], or in human cells using BAC arrays [25] did not provide the resolution to define replication domains and their boundaries. A tiling array study of human ENCODE regions covering 1% of the genome was also not able to precisely delineate replication domains, likely because they are typically larger than the 500 kb segments queried by this study [26,27]. Other high-resolution studies of specific chromosome segments in human and *Drosophila* [28-30] also did not delineate domain structure but noted that the relationship between replication timing and transcription was best described at the level of large multi-genic regions rather than individual genes [30].

We have recently mapped replication domain structure genome-wide in mouse embryonic stem cells (mESCs) and their differentiated counterparts [2]. We recognized the need for a comprehensive database to display these profiles and compare them between cell lines as well as to other chromosome-wide properties such as transcriptional activity or other epigenetic marks. While this can in principle be done on other public web browsers such as the UCSC Genome Brower http://genome.ucsc.edu, in many cases it is desirable to quickly visualize one's data from individual replicates or using unpublished data in a format that is not appropriate for public viewing. Such comprehensive public databases are becoming complex with the number of tissue and cell-type specific data sets that the reader must browse through. Generally speaking, they are tailored toward the display of static features of chromosomes, rather than dynamic cell-type specific features. Moreover, because of the increasing complexity of genome browsers, viewing large chromosomal regions necessary to visualize replication domains tends to be very slow (e.g. a 5-Mb chromosome segment takes several tens of seconds to display on the UCSC Genome Browser, but 2-3 seconds on ReplicationDomain).

For these reasons, we developed ReplicationDomain as a centralized repository that enables rapid comparative analysis of the genomic landscape for replication domain organization, with the potential to compare these properties to any other genome-wide chromosome data sets, such as those from ChIP-Chip or ChIP-Seq experiments. Our published microarray data sets are publicly available for any non-registered user to view and download. Furthermore, we have implemented a user registration system that allows registered users to upload their own data sets. Users have three options for data security, either to view their data sets privately without sharing (Über Private), share with other registered users (Private), or make their data sets publicly available on condition that they are published or "in press" in peer-reviewed journals (Public). In the future, we plan to implement additional data security mechanisms that will allow sharing of data sets only with designated registered users. The ability of registered users to upload data sets for private viewing facilitates confidential sharing of data prior to publication, or for quality control checks. ReplicationDomain uses an interactive interface designed to be intuitive for users familiar with the UCSC Genome Browser, with unique features optimized for rapid viewing of multi-megabase domain-wide chromosome properties and the option to jump to the same region of interest in the UCSC Genome Browser. Our recent demonstration of global re-organization of replication domains during ESC differentiation [2] predicts that replication profiling will provide a rapid and comprehensive means to evaluate cell-type specific features of global genome organization. ReplicationDomain will provide a valuable resource to consolidate replication profiles, making them more accessible to view and identify cell-type specific properties. We encourage users to begin uploading data sets and suggesting features that will improve this database.

Construction and content Current Contents

Currently, ReplicationDomain accommodates 10 microarray data sets publicly available without a login requirement. They consist of 7 replication profiles (genome-wide replication timing data) with probes spaced every 5.8 kb for three different mESC lines (D3, TT2 and 46C) either prior to differentiation (indicated ESC) or after differentiation to neural precursor cells (NPCs) either using defined medium in monolaver culture (indicated NPC/ASd6) or in conditioned medium as embryoid bodies (indicated NPC/EBM9), as well as one data set for iPS (induced pluripotent stem) cells. One replication timing profile has been performed with probes spaced every 100 bp along regions of chromosome 6 and 7 for undifferentiated mESC line D3. In addition, there are 2 transcription profiles for mESC line D3 and its NPC/EBM9 differentiated counterpart. The details of how these data sets were collected are explained under "Documentation" (Figure 1, button 1). Further details and conclusions derived from these data sets are presented in Hiratani et al [2]. Briefly, replication-timing data were obtained by hybridizing early and late replication intermediates simultaneously to Nimblegen oligonucleotide arrays containing evenly spaced probes (CGH arrays). The Data Display page shows the log₂ ratio of early to late replicating intermediates for each probe as a function of its position along the length of any given chromosome. Transcription data were obtained from Affymetrix arrays, using standard methods described under Documentation (Figure 1, button 1) or in Hiratani et al [2]. What is displayed is a bar, the height of which is the relative signal intensity, and the width of which delineates the start and stop sites of transcription. Gold-colored bars represent genes deemed sufficiently above background to be called "Present" (i.e. expressed) while garnet-colored bars represent genes considered to

be "Absent" (i.e. silent) or very lowly expressed. The web site is designed to anticipate any tab-delimited text data file as long as each record indicates the chromosome, start and end nucleotide coordinates, and a data value. Importantly, as long as the file structure follows this format, the data values can derive from any other genome-wide application such as ChIP-Chip, ChIP-Seq, and RNA-Seq.

System Architecture and Implementation

Software and Hardware

The ReplicationDomain software is an integrated system consisting of a relational database, PHP scripts, administrative interface, and an interactive web-based user interface (Figure 2). Upon upload as tab-delimited text, each data set is parsed and inserted into the database, which can then be queried for data in particular chromosomes, regions, or genes. PHP and HTML scripts are used to connect the user interface and uploaded data sets to the database, which is backed up twice daily to separate hard disk and weekly to tape. The database is hosted on MySQL Server 5.0.51, a multithreaded relational database management system. The Web based user interface, internal administrative interface, administrative scripts, and website statistics are programmed using Apache/2.2.9 (Free-BSD), mod_ssl 2.2.9, OpenSSL 0.9.8e DAV/2, and PHP 5.2.6 with Suhosin-Patch running on a Free BSD version 7.0 64-bit operating system. Graphical output is generated with ImageMagick 6.4.1.8 and GraphicsMagick 1.1.12, and access logs are recorded and analyzed using custom PHP scripts and Analog 6.0.1.

Utility and discussion

ReplicationDomain is designed to quickly and conveniently examine and compare properties of chromosomes that are important for their higher order structure and function, particularly as it relates to the organization of replication timing domains. It functions similar to other genome browsers, but has unique features that allow one to examine and compare microarray data for replication timing to steady state transcript levels (or any other genome-wide feature) with ease. Data can be downloaded as a tab-delimited text file or saved as a JPEG file for figure assembly. We also provide a link to the UCSC Genome Browser for quick access to additional data for the chromosomal region of interest. Members of the Gilbert lab utilize this site regularly as a tool to mine patterns in the genome indicative of functional changes in chromosome structure during differentiation (see "Downloading data for personal use"). Furthermore, we have implemented a user registration system that allows registered users to upload their own data sets, the details of which can be found under "Creating a ReplicationDomain account" and "Uploading your own data set" below. Regarding unpublished data sets, registered users can either view them privately without sharing or share with other regis-

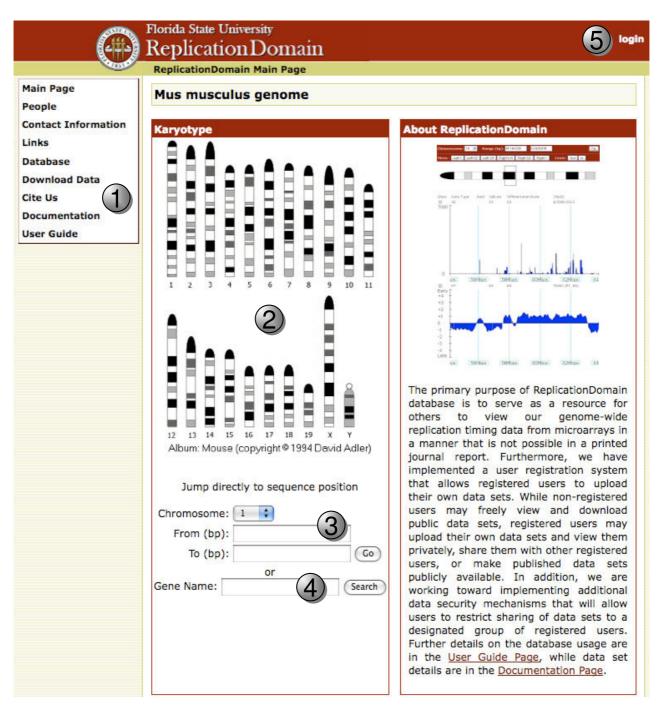


Figure I

ReplicationDomain Main Page. Buttons correspond to descriptions in the text.

tered users. For published or "in press" data sets, investigators have an option to release them for public access. We can accommodate data sets from any species with a physical genome map, but may need some time to assemble a new navigation path for each new species.

Database Interface

ReplicationDomain is freely accessible at <u>http://replica</u> <u>tiondomain.org</u>. It provides a novel tool to examine and compare microarray data for replication timing to other chromosomal properties. Different types of data sets can be aligned, by virtue of the dynamically generated graph-

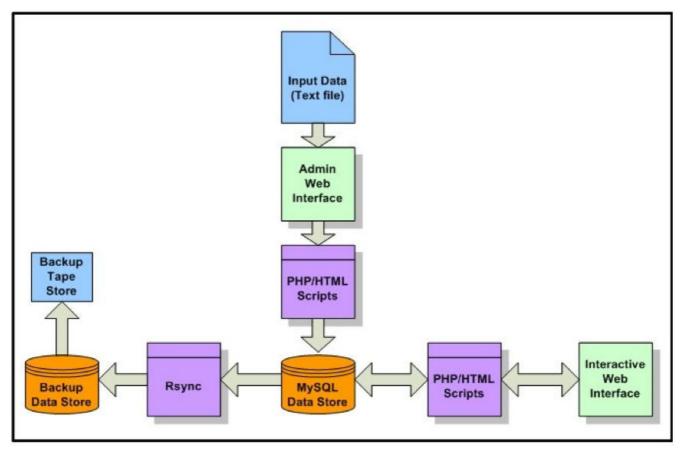


Figure 2

An Overview of ReplicationDomain Data Flow. Input data as text files is uploaded via administrative web interface, checked for errors, and inserted into MySQL database. Data is retrieved by interactive web interface for display. Data is archived via rsync to additional data store, data is written to tape from additional data store.

ical output. This unique feature allows investigators to compare replication timing and other chromosomal properties at any region of the genome under different developmental, genetic (e.g. gene knockout cell lines) and/or experimental conditions. Non-registered users may freely view and download public data sets. Registered users with a ReplicationDomain account may upload their own data sets and view them privately or share them with other registered users. Published or "in press" data sets can be added to the series of data sets public available. Further details of these procedures can be found below.

Getting to a region of interest

In the near future, species other than the mouse (*Mus musculus*) will be available and the user will need to first choose a species and be redirected to the data available for that species, similar to the UCSC Genome Browser. At present, only mouse data are uploaded so the viewer opens to the mouse page directly. From the Main Page (Figure 1), you can jump to the Data Display page (Figure

3) and view your favorite chromosomal region in one of three ways:

1) Click on the chromosome that contains your region on the ideogram (Figure 1, button 2). After being directed to the Data Display page of the chromosome, drag a rectangle around the region of interest on that chromosome (Figure 3, button 1).

2) Jump directly to a position by choosing the chromosome from the dropdown and typing in the desired base pair (bp) coordinates (Figure 1, button 3 and Figure 3, button 2).

3) Search for a Gene Name near the site of interest (Figure 1, button 4). You may be asked to choose from similarnamed genes, particularly if you have typed in a partial gene name. When you arrive to the chromosome position, your gene of interest will appear in the center of the image. You will likely need to zoom out to see the context of your gene. Your gene will remain at the center as you zoom.

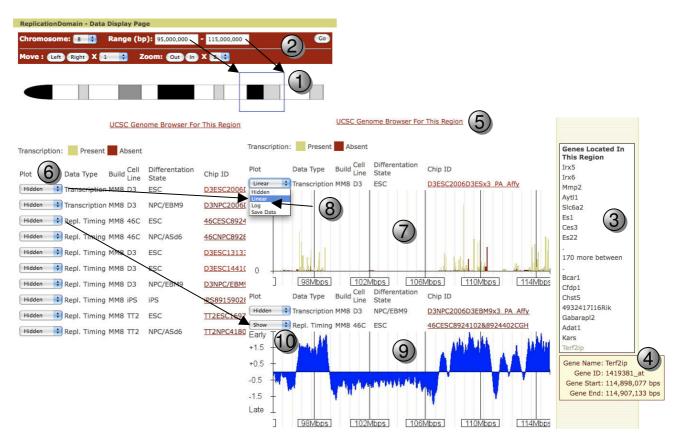


Figure 3 Data Display Page: An Exemplary View of A Region of Chromosome 8. Buttons correspond to descriptions in the text.

When you select a region of the chromosome, the Data Display page will display a floating gene name box (Figure 3, button 3) in the right column that contains the names and map positions of all genes within the selected region. If there are more than 16 genes, the first and last 8 genes are shown, with the number of genes between them indicated. Pointing your mouse cursor at a gene name will open a hover box that will provide its start and stop positions (Figure 3, button 4). Further information about this genomic region or the structure of the genes contained within it can be found by clicking the UCSC Genome Browser link (Figure 3, button 5).

Choosing data sets to view

Data sets on the Data Display page are chosen for viewing much as they are with the UCSC Genome Browser. All data sets are "Hidden" as a default, and can be viewed by using the dropdown menus (Figure 3, button 6). Relative transcription levels (Figure 3, button 7) can be plotted with either a linear or log scale (Figure 3, button 8) while replication-timing data is always plotted as a log_2 ratio (Figure 3, button 9) by choosing "Show" (Figure 3, button 10). When these options are selected, graphical display of

the data set will show up automatically. The y-axes are adjusted to the highest and lowest values in the entire data set (e.g. for transcription, the height of the y-axis provides an indication of the most highly transcribed gene in that data set). Details of each data set can be found by clicking on their Chip ID on the Data Display page (Figure 3, Chip ID column) or in the "Database" link in the main menu window (Figure 1, button 1). Definitions of the terms used to describe the data sets can be found under the "Documentation" link in the main menu (Figure 1, button 1).

Moving around

You can move to the left or the right, jump to a specific nucleotide position, or zoom in or out (2–8 fold) using the buttons at the top of the Data Display page (Figure 3, button 2), similar to what is done with the UCSC Genome Browser. ReplicationDomain also allows you to grab regions of the chromosome using your mouse cursor; simply drag a rectangle around the chromosomal region of interest either on the chromosome ideogram (Figure 3, button 1) or on individual data sets (Figure 3, buttons 7 and 9) and you will jump to that region. To orient your

position relative to genes in the area, you can: (1) use the floating gene name box in the right column (Figure 3, button 3); (2) use the link to the UCSC Genome Browser; or (3) download the corresponding "Transcription" data with all gene names and positions.

Downloading data for personal use

For further analysis (an example is illustrated in Figure 4), you can download the data sets shown on your screen (Figure 4A). When this option is chosen, you will be prompted to select a location on your desktop for downloading all the data points for that data set contained within the selected genomic region. The data will be

downloaded as a tab-delimited text file that will contain a description of the data at the top, followed by a table of data that you can move into any spreadsheet program of choice for further analysis (Figures 4B, 4C and 4D). You can download data representing a small region or up to an entire chromosome using this function. If you desire the complete data set for all chromosomes, that is provided separately through the "Download Data" link in the main menu window (Figure 1, button 1). Of course, you may also save any web page as a .pdf or take a snapshot using your own computer if you want to build quick figures for group presentations.

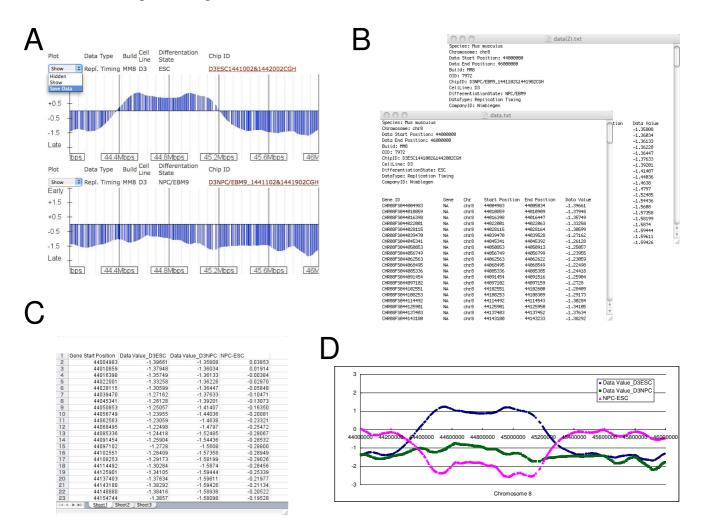


Figure 4

An Exemplary Analysis Performed with ReplicationDomain. Visual inspection with the aid of ReplicationDomain identifies a region where replication timing has changed during differentiation of ESCs to NPCs. Users can download the data shown on the screen by selecting the "Save Data" option of each data set (A). Downloaded data for each data set will be in a tab-delimited text file with a description at the top (B). Data from multiple differentiation states can be assembled and their differentials calculated in Microsoft Excel (C) and visualized graphically (D). The graph (D) shows a comparison of D3ESC (dark blue) and D3NPC/EBM9 (green) replication timing profiles and their differentials (NPC-ESC; magenta) on chr8:44,000,000– 46,000,000.

Comparing to other genome properties

At the top of the page is a link entitled "UCSC Genome Browser for this Region" (Figure 3, button 5). Clicking this button will open a new tab or browser page to the UCSC Genome Browser for the region viewed on ReplicationDomain. In addition, if you wish to upload your own data sets and compare to other data sets on ReplicationDomain, you will need to register yourself first (see "Creating a ReplicationDomain account" below for details). Registered users are entitled to upload their data sets as described below under "Uploading your own data set."

Creating a ReplicationDomain account

User registration is required for those who wish to upload their data sets or use the database interactively with other registered users. To create a ReplicationDomain account, please visit our ReplicationDomain Account Request page (a link is found on the User Guide page) and submit the Account Request Information form. A confirmation email will be sent to you with a ReplicationDomain username and a password, which will allow you to log in (Figure 1, button 5).

Uploading your own data set

To upload your data set, first generate a single, tab-delimited text file for each data set that includes the following 6 columns in this order as shown in Figures 5A and 5B: Gene/Probe ID (any unique identifier of genes or probes), Gene Name (type NA when unavailable), Chromosome, Chromosome Start, Chromosome End, and Data Value. You can also add a <u>Present-Absent</u> (P = 1, A = 0) column as a 7th column for transcription data sets. Then, log in with your ReplicationDomain account (Figure 1, button 5), go to the "Database" link from the main menu (Figure 1, button 1), click the "Upload data set" link and fill out the form (Figure 5C) with the requested attributes, select your text data file and hit "send." This should upload your data set, which will appear on the Data Display page for graphical display as well as on the list of data sets found under the "List my data sets" link on the Database page. All of the current 14 data set attributes are described under "Definitions of Data Entry Terms" below, which includes Data Security Level (description also on our Documentation page). Users can select Public, Private, or Über Private. Public data sets can be viewed by any user without a login requirement, but we require such data sets to be published or "in press" in peer-reviewed journals for database quality assurance purposes (i.e. a reference must be provided). Private data sets are viewable by all registered users with a ReplicationDomain account, while Über Private data sets are viewable only by the registered user who uploaded them. In the future, we anticipate implementing additional data security mechanisms that will allow sharing of data sets among designated registered users. Finally, existing data set attributes including Data Security Level can be edited later if necessary: click the "Database" link from the main menu, identify your data set from "List my data sets" link, hit "Edit," make necessary modifications, and hit "Update."

Definitions of Data Entry Terms

Data sets are defined by a combination of 14 entries, as described below (see also Figure 5C). Upon uploading data sets, registered users can either select terms from the dropdown list, or create a new term by filling in the blank.

Species

At present, we only have *Mus musculus*, but we intend to add *Homo sapiens* and *Drosophila melanogaster* soon. Contact us to create any new species page.

Company

Microarray product supplier name.

Chip ID

This is the unique identifier for each data set. While a "Chip ID" normally represents a single replicate experiment (e.g. one microarray hybridization), most data sets currently displayed on the Data Display page are averages of multiple replicates. Therefore, we have re-defined Chip ID as a string of characters combining individual "Chip ID" numbers and description of the data set identity. Chip ID is not useful except to communicate comments regarding a particular experiment. On the Data Display page, however, the Chip ID for each data set is set up as a link to the "Data Set Details" (Figure 3, Chip ID column; also accessible from the "Database" link on the main menu). The Chip ID is also useful for identifying data sets when downloading entire data sets through the "Download Data" link in the main menu.

Build

This indicates the version of the genomic sequence information that was used to assemble the microarray chip in the particular experiment (for example, mm7 or mm8 for the mouse). Builds change slightly as sequence information becomes updated, so the exact base pair position of any given DNA sequence will change as the sequence information becomes annotated. The build information indicated in each data set shows the build used for chromosomal coordinates of probes on the particular array type used.

Order ID

For Nimblegen data sets only.

Cell Line

This indicates the name of cell lines employed. At present, D3 [31], TT2 [32] AND 46C [33] are three different established mESCs. Α

CHR01FS003018759 NA chr1 3018759 3018888 CHR01FS003036253 NA chr1 3036253 303364 CHR01FS003036253 NA chr1 3041992 3042041 CHR01FS00305253 NA chr1 3041992 3042041 CHR01FS003052566 NA chr1 305566 3053655 CHR01FS003065156 NA chr1 3065156 3065265 CHR01FS003076536 NA chr1 3076536 307559 CHR01FS003076536 NA chr1 30976533 3093723 CHR01FS003105423 NA chr1 3105472 3093673 CHR01FS003105424 NA chr1 3105472 3105472 CHR01FS003116559 NA chr1 31105423 3105472 CHR01FS003124245 NA chr1 3142245 31322452 CHR01FS003156928 NA chr1 3162762 316631 CHR01FS003162762 NA chr1 3164281 3180487	Data
CHR01F3003036253 NA chr1 3036253 3036304 CHR01F5003041992 NA chr1 3041992 3042041 CHR01F5003041992 NA chr1 3045304 305365 CHR01F5003065366 NA chr1 3065156 3065205 CHR01F5003065156 NA chr1 3076536 3065205 CHR01F5003087094 NA chr1 3076536 3093722 CHR01F50030937053 NA chr1 3093673 3093722 CHR01F5003093673 NA chr1 305473 3093722 CHR01F5003105423 NA chr1 3105423 3105472 CHR01F5003116559 NA chr1 3116559 3116683 CHR01F500312452 NA chr1 3122452 3122504 CHR01F5003134254 NA chr1 3122452 3125094 CHR01F50031256928 NA chr1 3162762 316208 CHR01F5003127425 NA chr1 3162762 316208	-0.526375267100967
CHR01FS0030041992 NA chr1 3041992 3042041 LHR01FS003053606 NA chr1 3053605 3053665 3053665 LHR01FS003065156 NA chr1 3065165 3065565 3076598 LHR01FS0030676536 NA chr1 3065156 3076598 3076598 LHR01FS003087994 NA chr1 308794 3083972 LHR01FS003087994 NA chr1 3093722 3105423 LHR01FS003105423 NA chr1 3105423 3105472 LHR01FS003116559 NA chr1 3110559 3116688 LHR01FS003122452 NA chr1 3124252 3122504 LHR01FS003134224 NA chr1 3134224 3134594 LHR01FS0031355928 NA chr1 3165082 3156984 LHR01FS003162762 NA chr1 3162762 3162808 LHR01FS003162762 NA chr1 3162762 3162808 LHR01FS003162762 NA chr1 <td>-0.627372496310323</td>	-0.627372496310323
CHR01FS0030653606 NA chr1 3053606 3053655 CHR01FS003065156 NA chr1 3065156 3065255 CHR01FS003065156 NA chr1 3076536 3076536 CHR01FS003075536 NA chr1 3076536 3076536 CHR01FS003097537 NA chr1 3097653 3093773 CHR01FS003106423 NA chr1 3106542 3105472 CHR01FS0031065423 NA chr1 3116553 3116919 3110919 3110959 CHR01FS0031165559 NA chr1 3142542 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3122452 3162631 Chr1 3167262 3162631 Chr1 31679251 NA chr1 31678251 NA chr1 3168031 S16803 S168031 S168031 S168031 S168031 S168031	-0.7171644823638
CHR01FS003065156 NA chr1 3065156 3065205 CHR01FS003076536 NA chr1 3076536 3076536 3076536 3076536 3076536 3076536 3076536 30862051 CHR01FS003093673 NA chr1 3093673 3093722 3093722 CHR01FS003105423 NA chr1 3105423 3105472 CHR01FS003116559 NA chr1 3116553 3116608 CHR01FS003116559 NA chr1 3122452 3122452 CHR01FS0031156528 NA chr1 3122452 3125694 CHR01FS003156528 NA chr1 3146276 3156928 CHR01FS003155628 NA chr1 3156928 315694 CHR01FS003165628 NA chr1 3162762 3162762 CHR01FS003165628 NA chr1 3162762 3162762 CHR01FS003162762 NA chr1 3162762 3162802 CHR01FS003182551 NA chr1 3180251 318025	-0.743384080577782
CHR01FS003076536 NA chr1 3076536 3076598 CHR01FS003087994 NA chr1 308794 3088794 CHR01FS003087994 NA chr1 308794 3088795 CHR01FS003087994 NA chr1 308794 308825 CHR01FS003105423 NA chr1 3105423 3105472 CHR01FS003116559 NA chr1 3110559 3116688 CHR01FS003122452 NA chr1 3122452 312250 CHR01FS003134224 NA chr1 3134224 3134234 3134524 3134594 CHR01FS003155928 NA chr1 3165028 3156984 CHR01FS003162762 NA chr1 316762 316213 CHR01FS003162762 NA chr1 316762 31628251 CHR01FS003162762 NA chr1 3169263 3168036 CHR01FS003174487 NA chr1 3189251 3140835 CHR01FS003189508 NA chr1 3197183	-0.791542544482317
CHR01FS003087994 NA chr1 308794 3088051 CHR01FS003093673 NA chr1 3093673 3093723 CHR01FS003106423 NA chr1 3109542 309373 CHR01FS003106423 NA chr1 3105423 309372 CHR01FS003106525 NA chr1 3110559 3116968 CHR01FS003124252 NA chr1 312452 3122452 CHR01FS003134224 NA chr1 3134224 3134224 CHR01FS003134224 NA chr1 312452 3122452 CHR01FS003134224 NA chr1 3134224 3134273 CHR01FS003134224 NA chr1 3134224 314273 CHR01FS003134224 NA chr1 3162762 316281 CHR01FS003134224 NA chr1 3164251 3180305 CHR01FS0031274487 NA chr1 3180918 3186931 CHR01FS00312951 NA chr1 3189183 3186931	-0.832918085755317
CHR01FS003093673 NA chr1 3093673 3093722 CHR01FS003105423 NA chr1 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105423 3105424 NA chr1 3124523 3155943 3156943 3156943 3156943 3156943 3156943 3156943 3156943 3156943 3162762 NA chr1 3162762 3162762 3162763 3162763 3162763 3162616 NA chr1 3162763 3162613 3162613 317487 NA chr1 3180251 3180395 NA chr1 3180251 3180395 NA chr1 3180251 3180395 NA chr1 319155 3191583 3197183 3197183 3197183 3197183 3197183 <td>-0.8673162710887</td>	-0.8673162710887
CHR01FS003185423 NA chr1 3105423 3105423 CHR01FS003110519 NA chr1 3110919 3110960 CHR01FS003110559 NA chr1 3110559 3116580 CHR01FS003122452 NA chr1 3122452 312250 CHR01FS003124252 NA chr1 312452 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134234 3134244	-0.89555233366714
CHR01FS003110919 NA chr1 3110919 3110958 CHR01FS003122452 NA chr1 3116559 3116658 CHR01FS00312452 NA chr1 312455 312564 CHR01FS003134224 NA chr1 3134224 3134273 CHR01FS003134224 NA chr1 3156928 3134273 CHR01FS0031362762 NA chr1 3156928 3156928 CHR01FS0031362762 NA chr1 3156928 3162762 CHR01FS0031362762 NA chr1 3162762 3162762 CHR01FS00317487 NA chr1 3186926 3180305 CHR01FS003180251 NA chr1 3185980 3180305 CHR01FS003180251 NA chr1 3185980 316631 CHR01FS0031897138 NA chr1 3197187 3197187 CHR01FS003203197138 NA chr1 3203195 3208636 CHR01FS003214407 NA chr1 3203195 3208636	-0.907614473132328
CHR01FS003116559 NA chr1 3116559 3116608 CHR01FS003122452 NA chr1 3122452 3122545 DRR01FS003122452 NA chr1 3124224 3134224 DRR01FS003134224 NA chr1 3136252 3156948 DHR01FS003156928 NA chr1 3156928 3156948 DHR01FS003162762 NA chr1 3162762 3162762 DHR01FS003162762 NA chr1 3186276 3180305 DHR01FS003162762 NA chr1 3180262 3168305 DHR01FS003180251 NA chr1 3180305 3180305 DHR01FS003185508 NA chr1 3191525 319150 DHR01FS003197138 NA chr1 320197138 3197187 DHR01FS003208197138 NA chr1 32019713 32031252 DHR01FS003208197138 NA chr1 32019713 3208365 DHR01FS003208197138 NA chr1 3208365 3208656	-0.930614965070646
CHR01FS003122452 NA chr1 3122452 3122504 UHR01FS003134224 NA chr1 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134224 3134234 3162762 3162633 3168035 S168031 S168032 3180805 31808051 31808051 31808051 31808051 3197187 3197187 3197187 3197187 3197187 3197187 3197187 3197187 <t< td=""><td>-0.940590510275932</td></t<>	-0.940590510275932
CHR01FS003134224 NA chr1 3134224 3134224 CHR01FS003156928 NA chr1 3156928 3156928 CHR01FS003156928 NA chr1 3156928 3156928 CHR01FS003156926 NA chr1 3156928 3168916 CHR01FS003156762 NA chr1 3176479 3174549 CHR01FS003174487 NA chr1 3180251 3180305 CHR01FS0031880251 NA chr1 3180598 3186931 CHR01FS003185980 NA chr1 3180598 3186931 CHR01FS003191525 NA chr1 3197187 CHR01FS003203197138 NA chr1 3197137 CHR01FS003203197138 NA chr1 3283195 3208252 CHR01FS003203197138 NA chr1 3283195 3208252 CHR01FS003208195 NA chr1 3208636 3208656 CHR01FS003214407 NA chr1 3220863 3204646 CHR01FS0032214407 <t< td=""><td>-0.950389853285237</td></t<>	-0.950389853285237
CHR01FS003156928 NA chr1 3156928 3156984 CHR01FS003162762 NA chr1 3162762 3162761 CHR01FS003162762 NA chr1 3162762 3162761 CHR01FS003162762 NA chr1 3162762 3162761 CHR01FS003162762 NA chr1 3180261 3160305 CHR01FS003185580 NA chr1 3185980 3166031 CHR01FS003191525 NA chr1 3197187 3197187 CHR01FS003191525 NA chr1 3197138 3197187 CHR01FS003203195125 NA chr1 3203195 3208326 CHR01FS003203195135 NA chr1 3201453 3208456 CHR01FS003204636 NA chr1 3204636 3208656 CHR01FS003214407 NA chr1 3220863 3214467 CHR01FS003220464 NA chr1 3220804 32240653	-0.960220616553999
CHR01FS003162762 NA chr1 3162762 3162816 CHR01FS003174487 NA chr1 3174487 3174487 PHR01FS003174487 NA chr1 3174487 3174487 PHR01FS003174487 NA chr1 3174487 3104305 PHR01FS003184251 NA chr1 3185980 3186031 CHR01FS003191525 NA chr1 3191525 3191580 PHR01FS003203195125 NA chr1 3191525 3191580 CHR01FS003203195125 NA chr1 3191525 3191580 CHR01FS003203195135 NA chr1 3201636 3208636 CHR01FS0032031955 NA chr1 3208636 3208636 CHR01FS003203195 NA chr1 3208636 3208636 CHR01FS003203195 NA chr1 3208636 3208636 CHR01FS0032208047 NA chr1 3208636 3208696 CHR01FS003220804 NA chr1 3220804 32208636	-0.978883518109678
CHR01FS003174487 NA chr1 3174487 3174549 CHR01FS003180251 NA chr1 3180251 3180365 CHR01FS003185980 NA chr1 3180263 3186985 CHR01FS003185980 NA chr1 3180263 3186985 CHR01FS003195980 NA chr1 3191525 3191586 CHR01FS003203197138 NA chr1 3197138 3197187 CHR01FS0032031955 NA chr1 3203195 3208252 CHR01FS0032080636 NA chr1 3208045 3208056 CHR01FS003214407 NA chr1 3224407 3214407 CHR01FS0032208064 NA chr1 3220804 32208053	-1.01299052320285
CHR01FS003180251 NA chr1 3180251 3180305 CHR01FS003185980 NA chr1 3185980 3186931 CHR01FS003185980 NA chr1 3185980 3186931 CHR01FS003191525 NA chr1 3197152 319153 CHR01FS003197138 NA chr1 3197133 3197187 CHR01FS003203195 NA chr1 3203195 32083252 CHR01FS0032080536 NA chr1 3208636 3208656 CHR01FS003214407 NA chr1 3214407 3214407 CHR01FS003220004 NA chr1 3220004 3220053	-1.01967905722572
CHR01FS003185980 NA chr1 3185980 3186031 LHR01FS003191525 NA chr1 3191525 319152 LR01FS003191525 NA chr1 3191525 319150 LR01FS003197138 NA chr1 3197138 3197137 CHR01FS003203195 NA chr1 3208195 3208252 LHR01FS003208636 NA chr1 3208636 3208656 LHR01FS003214407 NA chr1 3208636 32086964 LHR01FS003220804 NA chr1 3220804 3220804	-1.02415565436567
CHR01FS003191525 NA chr1 3191525 3191580 CHR01FS003197138 NA chr1 3197138 3197187 CHR01FS003203195 NA chr1 3203195 3203195 3203195 3208252 CHR01FS003208195 NA chr1 3208195 3208265 Chr1 3208195 3208266 CHR01FS003214407 NA chr1 3214407 3214467 3214467 3214467 CHR01FS003220004 NA chr1 3220804 32208053 3208053	-1.02405856899875
CHR01FS003197138 NA chr1 3197138 3197187 CHR01FS003203195 NA chr1 3203195 3203252 CHR01FS0032080636 NA chr1 3208636 3208053 CHR01FS003208124807 NA chr1 3214407 3214407 CHR01FS0032208064 NA chr1 3220804 3220804	-1.02382167238355
CHR01FS003203195 NA chr1 3203195 3203252 CHR01FS003208636 NA chr1 3208636 3208636 CHR01FS003208636 NA chr1 3208636 3208636 CHR01FS003208637 NA chr1 3214407 3214467 CHR01FS003220004 NA chr1 3220823 3220804	-1.02440692577355
CHR01FS003208636 NA chr1 3208636 3208696 CHR01FS003214407 NA chr1 3214407 3214464 CHR01FS003220004 NA chr1 3220004 3220053	-1.02675658239546
CHR01FS003214407 NA chr1 3214407 3214464 CHR01FS003220004 NA chr1 3220004 3220053	-1.03235820306277
CHR01FS003220004 NA chr1 3220004 3220053	-1.04178660772096
	-1.05680686794487
	-1.07525426015840
	-1.14644803422481
	-1.16968953591672
	-1.21193000075663 -1.22674361368429

D3ESC_Transcription.txt

Start End Data 4334225 4350473 10.6

4763291 4775791 1933.2

4797973 4836815 1627.2

4847894 4887983 2208.4

4900049 5009460 18.9

5073253 5152630 452.3

5578573 5592947 13.7

6720131 6851021 9.4

8352594 9284944 5.3

6204742 6265656 126.6

7079000 7163563 296.2

9530655 9532665 983.8

9533293 9563182 20.2

4481008 4486494 25

Present-Absent

1

Chromosome End Position (Integer number) Chromosome End Position (Integer number) Data (Floating point number) Present or Absent (This is an optional column. "1" = present, "0" = absent) Contact us to request additional species Species: Company: Previously used values 🗦 Chip ID: Previously used values Previously used values 🗦 Build: Order ID: Previously used values Cell Line: Previously used values Differentiat Previously used values State: Array Desig Name: Previously used values \$ Replication Timing - Contact us to request additional data types Data Type: Reference Comments: Yes 🕘 No 💿 Present or Absent column included?: Public O Private O Uber Private 💿 Data Security Level: • Public: Dataset is publicly viewable Private: Only registered users can view the dataset.
Uber Private: Only you can view the dataset.
A reference must be included to make Public datasets publicly viewable. Data starts on line: 1 🛟 Upload data file: 選択... Send File

ReplicationDomain - Database - Data Upload Columns should appear in the following order in a tab delimited text file

ID (Probe or datapoint identifier)
 Gene Name (Use "NA" if unknown)
 Chromosome (ex, "chr1")

Figure 5

B

000

1424964_at

1429177_x_at 1425189_a_at

1448244_at

1419258_at

1443694_at

1415826_at

1451813_at

1418968_at

1436793 at

1435634_at

1416998_at

1424392_at

1432287_a_at

Gene

Rp1h

Sox17

Mrpl15

Lypla1 chr1

Tcea1

Rgs20

Oprk1

St18

Rb1cc1 chr1

Pemtd1 chr1

Sntg1

Adhfe1 chr1

Rrs1

Atp6v1h chr1

chr

chr1

chr1

chr1

chr1

chr1

chr1

chr1

chr1

chr1

AffyID

Uploading Your Own Data Set. (A, B) Data file structures. Exemplary tab-delimited text data files for replication timing (A) and transcription (B). (C) A screenshot of the "Data Upload" page.

Ŧ

Differentiation State

This is the tissue or tissue type represented by the cell line used and grown under the indicated conditions. Currently, there are four such differentiation states:

ESC: Undifferentiated ES cells

NPC/EBM9: The 9th day of differentiation following an established neural differentiation protocol that differentiates ESCs via embryoid bodies to Sox1 positive NPCs in conditioned medium [34]. NPC/ASd6: The 6th day of differentiation following an established neural differentiation protocol that differentiates ESCs to Sox1 positive NPCs in monolayer cultures using defined medium [33].

iPS: "Induced Pluripotent Stem Cells" re-programmed from tail-tip fibroblasts derived from a 129xBL-6 hybrid strain of mice to the pluripotent state as described [35].

Array Design Name

Microarray supplier and catalog number.

Data Type

Indicates the property being measured in the indicated experiment. At present, replication timing and transcription data are shown. In the future, data for other genomewide properties of chromosomes may be displayed. Contact us to add new types of data such as ChIP-Chip or ChIP-Seq.

Reference

Data sets to be displayed publicly must include a reference.

Comments

We provide detailed microarray design information here but any additional comments can be added.

Present or Absent Column

For uploading transcription data sets that contain presentabsent calls, specify here.

Data Security Level

Users can select Public, Private, or Über Private. Users can make their published or "in press" data sets publicly available by selecting "Public" and providing a reference under the entry term, *Reference*. Private data sets are viewable by all registered users with a ReplicationDomain account, while Über Private data sets are viewable only by the user who uploaded the data set.

Data Starts on Line

Usually starts on line 2, with line 1 being the column names.

Conclusion

ReplicationDomain provides a user-friendly platform to view replication timing data from any organism and compare that data to other properties in a manner that is optimized for rapid viewing of multi-megabase segments of chromosomes. In addition to providing a consolidated and devoted database, ReplicationDomain provides the opportunity for researchers to share and analyze preliminary data sets with colleagues prior to providing public access. Although not as comprehensive as other databases, ReplicationDomain allows rapid linkage to the UCSC Genome Browser for cross-referencing to other databases. At present, the site contains only data sets collected in the Gilbert laboratory. In the future, we expect to serve as curators of a substantial database, and we expect ReplicationDomain to contain data from other laboratories and for other species, as well as other specific properties that are relevant to higher order domain structure of chromosomes. We invite others to use the web site and to create an account and upload their own data sets so that ReplicationDomain can be used to advance our understanding of the functional significance of a dynamic, developmentally regulated replication-timing program.

Availability and requirements

ReplicationDomain is available at <u>http://www.replica</u> <u>tiondomain.org</u> for use by academic or non-academic users without restriction or charge.

Authors' contributions

NW and AS are the web site designers, NW designed and programmed interactive web interface and updates/maintains administrative interface, IH generated the data and consulted in all aspects of the design, TR and TY assisted with implementation of the design and communications between biologists and computer scientists, AS designed hardware and software architecture, designed database structure, designed and programmed administrative interface and scripts, DMG conceived of the project and coordinated all aspects of its development. IH, AS, and DMG wrote the manuscript and all authors have read and approved the manuscript.

Acknowledgements

We thank S. Thompson (FSU) for technical lectures and suggestions in initial sessions. D. Adler (ZymoGenetics, University of Washington) for mouse chromosome ideograms. This work was supported by NIH award GM83337 to DMG and postdoctoral fellowships from the Rett Syndrome Research Foundation to IH and the Leukemia & Lymphoma Society to TY.

References

- Gilbert DM, Gasser SM: Nuclear Structure and DNA Replication. In DNA replication and human disease Edited by: DePamphilis ML. Cold Spring Harbor, New York: Cold Spring Harbor Press; 2006.
- Hiratani I, Ryba T, Itoh M, Yokochi T, Schwaiger M, Chang CW, Lyou Y, Townes TM, Schubeler D, Gilbert DM: Global reorganization of replication domains during embryonic stem cell differentiation. PLoS Biol 2008, 6(10):e245.
- Gilbert DM: Replication timing and transcriptional control: beyond cause and effect. Curr Opin Cell Biol 2002, 14:377-383.
- Goren A, Cedar H: Replicating by the clock. Nat Rev Mol Cell Biol 2003, 4(1):25-32.
- Schwaiger M, Schubeler D: A question of timing: emerging links between transcription and replication. Curr Opin Genet Dev 2006, 16(2):177-183.
- Loupart M, Krause S, Heck MS: Aberrant replication timing induces defective chromosome condensation in Drosophila ORC2 mutants. Curr Biol 2000, 10(24):1547-1556.
- Pflumm MF, Botchan MR: Orc mutants arrest in metaphase with abnormally condensed chromosomes. Development 2001, 128(9):1697-1707.
- Amiel A, Elis A, Sherker S, Gaber E, Manor Y, Fejgin MD: The influence of cytogenetic aberrations on gene replication in chronic lymphocytic leukemia patients. Cancer Genet Cytogenet 2001, 125(2):81-86.
- Sun Y, Wyatt RT, Bigley A, Krontiris TG: Expression and replication timing patterns of wildtype and translocated BCL2 genes. Genomics 2001, 73(2):161-170.
 Smith L, Plug A, Thayer M: Delayed replication timing leads to
- Smith L, Plug A, Thayer M: Delayed replication timing leads to delayed mitotic chromosome condensation and chromosomal instability of chromosome translocations. Proc Natl Acad Sci USA 2001, 98(23):13300-13305.
- Amiel A, Elis A, Maimon O, Ellis M, Herishano Y, Gaber E, Fejgin MD, Lishner M: Replication status in leukocytes of treated and untreated patients with polycythemia vera and essential thrombocytosis. Cancer Genet Cytogenet 2002, 133(1):34-38.
- Korenstein-Ilan A, Amiel A, Lalezari S, Lishner M, Avivi L: Allele-specific replication associated with aneuploidy in blood cells of patients with hematologic malignancies. Cancer Genet Cytogenet 2002, 139(2):97-103.
- Breger KS, Smith L, Thayer MJ: Engineering translocations with delayed replication: evidence for cis control of chromosome replication timing. *Hum Mol Genet* 2005, 14(19):2813-2827.
- Chang BH, Smith L, Huang J, Thayer M: Chromosomes with delayed replication timing lead to checkpoint activation, delayed recruitment of Aurora B and chromosome instability. Oncogene 2007, 26(13):1852-1861.

- Barbosa AC, Otto PA, Vianna-Morgante AM: Replication timing of homologous alpha-satellite DNA in Roberts syndrome. Chromosome Res 2000, 8(7):645-650.
- Hansen RS, Stoger R, Wijmenga C, Stanek AM, Canfield TK, Luo P, Matarazzo MR, D'Esposito M, Feil R, Gimelli G, et al.: Escape from gene silencing in ICF syndrome: evidence for advanced replication time as a major determinant. Hum Mol Genet 2000, 9(18):2575-2587.
- Hellman A, Rahat A, Scherer SW, Darvasi A, Tsui LC, Kerem B: Replication delay along FRA7H, a common fragile site on human chromosome 7, leads to chromosomal instability. *Mol Cell Biol* 2000, 20(12):4420-4427.
- Reish O, Gaí R, Gaber E, Sher C, Bistritzer T, Amiel A: Asynchronous replication of biallelically expressed loci: a new phenomenon in Turner syndrome. Genet Med 2002, 4(6):439-443.
- Norio P, Kosiyatrakul S, Yang Q, Guan Z, Brown NM, Thomas S, Riblet R, Schildkraut CL: Progressive activation of DNA replication initiation in large domains of the immunoglobulin heavy chain locus during B cell development. *Mol Cell* 2005, 20(4):575-587.
- Farkash-Amar S, Lipson D, Polten A, Goren A, Helmstetter C, Yakhini Z, Simon I: Global organization of replication time zones of the mouse genome. *Genome Res* 2008, 18(10):1562-1570.
- 21. Watanabe Y, Fujiyama A, Ichiba Y, Hattori M, Yada T, Sakaki Y, Ikemura T: Chromosome-wide assessment of replication timing for human chromosomes II q and 21 q: disease-related genes in timing-switch regions. *Hum Mol Genet* 2002, II(1):13-21.
- Watanabe Y, Ikemura T, Sugimura H: Amplicons on human chromosome IIq are located in the early/late-switch regions of replication timing. *Genomics* 2004, 84(5):796-805.
- MacAlpine DM, Bell SP: A genomic view of eukaryotic DNA replication. Chromosome Res 2005, 13(3):309-326.
- Schubeler D, Scalzo D, Kooperberg C, Van Steensel B, Delrow J, Groudine M: Genome-wide DNA replication profile for Drosophila melanogaster: a link between transcription and replication timing. Nat Genet 2002, 32(3):438-442.
- Woodfine K, Fiegler H, Beare DM, Collins JE, McCann OT, Young BD, Debernardi S, Mott R, Dunham I, Carter NP: Replication timing of the human genome. Hum Mol Genet 2004, 13(2):191-202.
- Karnani N, Taylor C, Malhotra A, Dutta A: Pan-S replication patterns and chromosomal domains defined by genome-tiling arrays of ENCODE genomic areas. Genome Res 2007, 17(6):865-876.
- Birney E, Stamatoyannopoulos JA, Dutta A, Guigo R, Gingeras TR, Margulies EH, Weng Z, Snyder M, Dermitzakis ET, Thurman RE, et al.: Identification and analysis of functional elements in 1% of the human genome by the ENCODE pilot project. Nature 2007, 447(7146):799-816.
- White EJ, Emanuelsson O, Scalzo D, Royce T, Kosak S, Oakeley EJ, Weissman S, Gerstein M, Groudine M, Snyder M, et al.: DNA replication-timing analysis of human chromosome 22 at high resolution and different developmental states. Proc Natl Acad Sci USA 2004, 101(51):17771-17776.
- Woodfine K, Bearé DM, Ichimura K, Debernardi S, Mungall AJ, Fiegler H, Collins VP, Carter NP: Replication Timing of Human Chromosome 6. *Cell Cycle* 2005, 4(1):172-176.
 MacAlpine DM, Rodriguez HK, Bell SP: Coordination of replica-
- MacAlpine DM, Rodriguez HK, Bell SP: Coordination of replication and transcription along a Drosophila chromosome. Genes Dev 2004, 18(24):3094-3105.
- Doetschman TC, Eistetter H, Katz M, Schmidt W, Kemler R: The in vitro development of blastocyst-derived embryonic stem cell lines: formation of visceral yolk sac, blood islands and myocardium. J Embryol Exp Morphol 1985, 87:27-45.
- Yagi T, Tokunaga T, Furuta Y, Nada S, Yoshida M, Tsukada T, Saga Y, Takeda N, Ikawa Y, Aizawa S: A novel ES cell line, TT2, with high germline-differentiating potency. Anal Biochem 1993, 214(1):70-76.
- Ying QL, Stavridis M, Griffiths D, Li M, Smith A: Conversion of embryonic stem cells into neuroectodermal precursors in adherent monoculture. Nat Biotechnol 2003, 21(2):183-186.
- Rathjen J, Haines BP, Hudson KM, Nesci A, Dunn S, Rathjen PD: Directed differentiation of pluripotent cells to neural lineages: homogeneous formation and differentiation of a neurectoderm population. Development 2002, 129(11):2649-2661.

